


# Safety Easy Temperature checker -Onetsu no Mihariban-

Safety Easy Temperature

Checker 37.5

## Onetsu no Mihariban

Different Colors appear based on your body temperature.

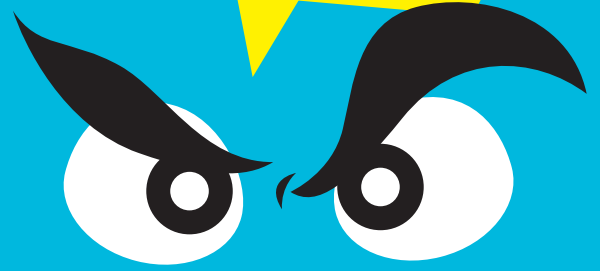

This sticker reminds you of your **body temperature** to **release the fear of high-fever** during Corona period.


Managers can monitor health condition of their staffs and **give good impression to your clients** because of **self-monitoring management** with the result of **higher CS**.

### Easy to use!

You can see the current fever condition based on colors



**Low Fever**  
(less than 33 degree)

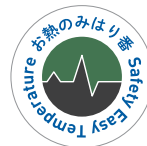

**Normal Fever**  
(between 35 and 37.5)



**High Fever**  
(more than 37.5 degree)


How to use ; Stick that on your body and wait & see!

1



Peel a sticker off the sheet

2

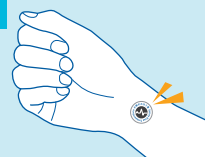

Stick it on your hand and wait a bit

3



Confirm the color based on body temperature

# Safety Easy Temperature checker -Onetsu no Mihariban-

## Safety Instruction

### Before you use

- Check that the sticker is not torn or damaged.
- Remove body warmers and coolants that affect body temperature measurement from your body before using the sticker.

### How to use

- Attach it to a place where body temperature can be easily measured (wrist, under ear, forehead, neck, etc.). It is recommended to attach it to the neck under the ear where the temperature change is most sensitive.
- Wipe off sebum, sweat, and makeup on the part where the sticker is attached.
- From a hygienic point of view, we recommend using within 24 hours.

### Storage method

- Store in a closed container indoors (30° C or less), avoiding high temperature and humidity and direct sunlight.

### Cautions

- The sealed thermometer roughly determines the body temperature.
- There may be an error in the color change due to the outside temperature and environment.
  - **Do not put it on scratches.**
  - **Although it is harmless to the human body, do not put it in your mouth.**
  - **Should you experience any skin problems (rash, redness, itchiness, etc.) due to use, discontinue use immediately and consult a doctor.**
  - **When using for small children, please be careful when using and storing as there is a risk of accidental ingestion.**

### Product info

Use adhesives that are safety certified according to RoHS standards.

1 Package : Contains 12 stickers

<https://www.set-corp.co.jp/mihariban/sticker.html>

Front



Back

